

- These two items can be used for various things, of which the subject ID is the most obvious and most impacting. So, when importing subjects from a previous study (either from Viedoc or an external system), it is important that these are set properly, especially if they have been used as part of the subject ID and you want to keep the same subject ID.
- If those extensions are present, their content is imported.

Note! If the sequence numbers exist in the ODM file, but they are also mapped during the import process (at Step3/5), then the mapping takes precedence.

- If those extensions are not present in the imported file, it is possible to specify them during the import process (Step 3/5).
- When importing to an existing study with existing data, matching of existing subjects is performed as follows:
 - If the sequence numbers are provided, these are used to perform the matching.

Note! If the sequence numbers exist in the ODM file, but they are also mapped during the import process (at Step3/5), then the mapping takes precedence.
 - If the sequence numbers are provided, but no match is found then the `SubjectKey` is used.
- If the sequence numbers are not present, then one of these needs to be mapped during the import process (Step 3/5).


5 Workflow - Import a new study from ODM

- Before you start importing the ODM file, you have to make sure that you already have a study in Viedoc that has a study design that matches the data structure in the ODM file to be imported. The metadata version(s) in the ODM file to be imported must contain all the events, forms, itemgroups, items and code list values that are referenced by ClinicalData. The import process performs the matching only by using OIDs.

In case you don't have such a study, this can be done as described below:

- Create the study in Viedoc Admin and make sure that you invite a user to the study that would have access to Viedoc Designer to perform the design import mentioned below (for instructions see [Add new study](#)). **(A)**
- In Viedoc Designer, import the design from the ODM file to the study you have just created previously in Viedoc Admin (for instructions see [Initiating a design](#)) **(B)**
- **Note!** This step must be performed only if the ODM file to be imported does not originate from Viedoc.

In Viedoc Designer, for the newly imported study design, go to **Study Settings** and configure the **Subject ID Generation Settings** (for instructions see [Study settings](#)). This will impact the selection you would have to make later on in *Import ODM File - Step 3/5* of the ODM import described in the later sections. **(C)**

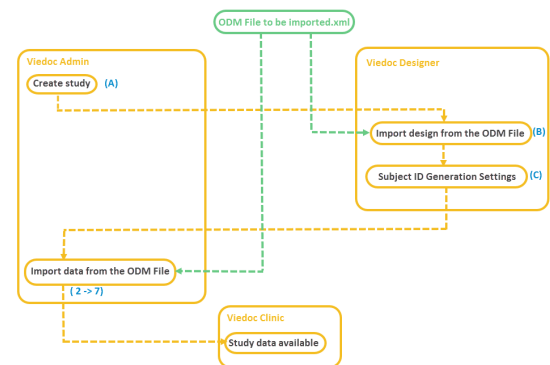

6 Initiate the import

- Select your study where the data will be imported
- Go to **Study Settings**
- Select the **Import ODM File** tab

7 Import ODM File - Step 1/5

- Click the **Upload a file** button and browse for your ODM file to be imported.
- The file name and size will appear right under the **Upload a file** button. You might get some error messages in red saying that the file cannot be uploaded due to some missing content (according to the CDISC ODM standard). In this case you have to go back to your ODM file, fix the error and upload the file again.
- Click **Upload and continue**. This will get you to the next step 2/5 of the import process.

Study settings
Here you can set settings for study.

Settings Date & time format Medical Coding **Import ODM File** API configuration

Import ODM File Step 1/5

Upload a file

DEMO_123_20170503_090338.xml 0.3MB

Import to demo

Upload and continue

History

File Name	Date and Time	User Name	Status	Log
-----------	---------------	-----------	--------	-----

8 Import ODM File - Step 2/5

- Choose from the *Study design version* dropdown list, the version of your study design in Viedoc the data will be imported to. The design has to match exactly your ODM data to be imported.
- In **Study site mappings**, the two columns on the left represent all the sites found in the ODM file you have uploaded. On the right side column are the sites currently available in the study you have selected to import to.

The system will perform an automatic mapping based on site name (not case sensitive). If the *code* extension is present (if the ODM file originates from Viedoc), this is mapped as well. If this is empty, only the name is mapped.

If no match is found it will default map to "Create new site". You can also manually perform the mapping by selecting in the drop-down list. If a match is found but you anyway choose *Create new site*, a duplicate site will be created, which is not a recommended approach.

Note! Make sure that every *Location* in the ODM file to be imported has at least one *MetaDataVersionRef* defined, otherwise no designed will be assigned to the respective site.

- Enter an email address to be used when the import needs to create audit-records, in the *Email address to be used for invented audit records*.
- Click **Continue**. This will get you to the next step 3/5 of the import process.

Study settings
Here you can set settings for study.

Settings Date & time format Medical Coding **Import ODM File** API configuration

Import ODM File Step 2/5

DEMO_123_20170503_090338.xml - Source system: VIEDOC 4.34.6310.23264 - Subjects: 4

Metadata version to study design version mappings

Metadata version OID (from xml) Study design version

1.0 Demo study 2014 2.0

Study site mappings

LOC:2815 Uppsala Uppsala (02)

E-mail address to be used for invented audit records

audit@viedoc.com

Continue

9 Import ODM File - Step 3/5

- Under **Study event dates** choose what dates you want to be matched to you events. If no selection is made and a form and item combination within the event called "\$EVENT.EventDate" is found (the Viedoc way to store event dates), this will be used. If the ODM file originates from Viedoc, you will find this form/item combination in the drop-down-list.
- From the **Populate** dropdown lists, you can set how the subjects will be ordered in Viedoc, according to the order they have in the ODM file to be imported. It is important to do this correctly, because the settings here determine how the subjects will be ordered and matched in both cases when no data exist in Viedoc and when the study already has data in Viedoc and you import some more (in this case it is important that the subjects are matched correctly so that no duplicates are created).

Study settings
Here you can set settings for study.

Settings Date & time format Medical Coding **Import ODM File** API configuration

Import ODM File Step 3/5

DEMO_123_20170503_090338.xml - Source system: VIEDOC 4.34.6310.23264 - Subjects: 4

Study event dates

Add subject SEVENT/EventDate

Visit 1 SEVENT/EventDate

Visit 2 ---

Home adm. ---

Visit 3 VS/VSDT

Unscheduled EC/ECDT

Medical / Surgical History STAT/STATDT

Prior and Concomitant Medications STAT/STATDDO

Adverse Events STAT/STATWDDT

SEVENT/EventProposedDate

SEVENT/EventPlannedDate

SEVENT/EventWindowStartDate

SEVENT/EventWindowEndDate

SEVENT/EventDate

Populate Populate

--- First ---

Back Continue

When a subject is added in Viedoc, there are two system items created for it:

- *SiteSubjectSeqNo*

- *StudySubjectSeqNo*

These items can be used, for example, when setting the Subject ID (in **Subject ID Generation Settings** mentioned in section 1. *ODM import overview*).

The settings to be performed here depend on if the ODM file to be imported has the *SiteSubjectSeqNo* and/or *StudySubjectSeqNo* extensions (ODM file originates from Viedoc) or not.

- If the ODM file has the *SiteSubjectSeqNo* and/or *StudySubjectSeqNo* extensions (ODM file originates from Viedoc), then no settings are required at this step, as they will be imported. You can continue to Step 4/5.
 - In case the ODM file does not have the *SiteSubjectSeqNo* and/or *StudySubjectSeqNo* extensions, this step depends on the settings made in Viedoc Designer regarding the **Subject ID Generation Settings** (described in the first section 1. *ODM import overview*).
 - If you have used one of the *SiteSubjectSeqNo* or *StudySubjectSeqNo*, then you have to choose this one from the first **Populate** dropdown list, and then choose from which form in your study this variable will be picked from in the **(form/item)** dropdown list. In the second **Populate** dropdown list you can choose if the first or last occurrence in the data set should be used, if the variable you point out occurs multiple times for a single subject (default is *First*).
 - If you have used a different one than the *SiteSubjectSeqNo* and *StudySubjectSeqNo*, then you should not make any selection in the **Populate** fields.
- **Note!** Once you have selected an option from the dropdown list, it is not possible to go back to no selection. It is only possible to select another option from the list.
 - Click **Continue**. This will get you to the next step 4/5 of the import process.

10 Import ODM File - Step 4/5

- **Select events to be excluded** by clicking on the field and choosing the events from your ODM file that you don't want to be included in the imported study. You can select one event at a time and repeat the operation for all the events that you don't want to include in the import.
- **Select forms to be excluded** by clicking on the field and choosing the forms from your ODM file that you don't want to be included in the imported study. You can select one form at a time and repeat the operation for all the forms that you don't want to include in the import.

Note! When importing an ODM file that was exported from Viedoc 4, you must exclude the *\$EVENT* form.

- Click **Continue**. This will get you to the last step 5/5 of the import process.

11 Import ODM File - Step 5/5

- A list of the users (from the imported ODM file) identified by email address and full name is displayed.
- Confirm the list of the users to be added to your study by entering your password and clicking **Import**.
- **Note!** The users are not active at this moment. You have to send invitations to these users using these email addresses in order for them to get login access to Viedoc (for instructions see [User Management](#)).

12 After the import

- Once you have imported the file and invited the users, your study is now in Viedoc, accessible for the users you have invited according to the respective roles.

- The PDFs with the forms history are not available right after the import, but once you have performed and export to PDF (for instructions see [Data export](#)), those will be generated and available in Viedoc.